



Disposable Laparoscopic Specimen Bags

The self-opening, high-strength bag ensures safe specimen capture with easy, one-handed operation.




UNIVERSAL PLUS
SUCTION/IRRIGATION



DETACHATIP® III
SYSTEM



VCARE®
RETRACTOR/ELEVATOR

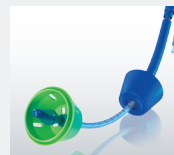

To learn more about this and other innovative products, call 800-448-6506 or visit ConMed.com.



525 French Road, Utica, NY 13502

DISPOSABLE
LAPAROSCOPIC
SPECIMEN
BAGS

The ConMed EndoSurgery Specimen Bag System is a sterile, disposable device designed for use as a receptacle for the collection and extraction of tissue, organs, and calculi during laparoscopic surgical procedures. The Specimen Bag System consists of a flexible memory wire specimen bag and introducer. The Specimen Bag System can be used through a 10, 11, or 12 mm cannula.

Unique Design

The unique patented design ensures ease of use, and can be deployed and manipulated with one hand. The self-opening memory wire, high-strength large bag offers easy and safe specimen capture while resisting tearing, reducing the potential for cell migration.

Easy to Use

Indicators on the thumb ring position the opening (top) of the bag for correct placement. Bag can be detached within abdominal cavity for later retrieval.


Large bag volume to capture large specimens

Minimum insertion force enables smooth deployment

Ergonomic design provides easy one-handed operation

Insert indicator ensures bag will be in proper position when deployed

String with slip knot facilitates closure of the specimen bag

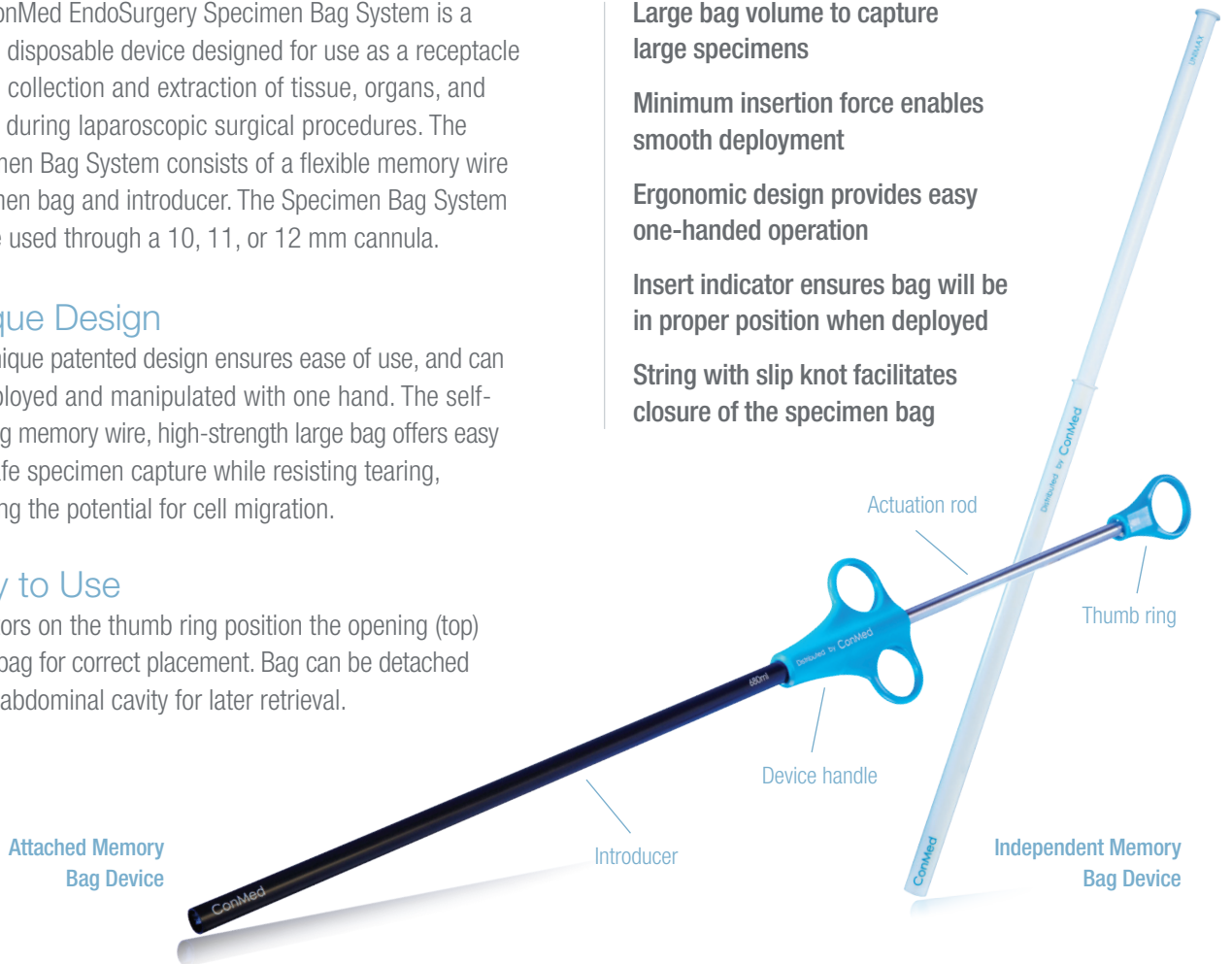

Disposable pouch has 510K clearance and is CE Marked.

Ordering Information

Description	Diameter	Length	Volume	Cannula	Packaging	Catalog Number
Independent Memory Wire Bag	3in (7.62cm)	6in (15.24cm)	250ml	10-12mm	1/pkg, 10/cs	SB836
Independent Memory Wire Bag	5in (12.7cm)	7in (17.78cm)	500ml	10-12mm	1/pkg, 10/cs	SB857
Attached Memory Wire Bag	3in (7.62cm)	6in (15.24cm)	250ml	10-12mm	1/pkg, 5/cs	SB936
Attached Memory Wire Bag	5in (12.7cm)	7in (17.78cm)	500ml	10-12mm	1/pkg, 5/cs	SB957
Attached Memory Wire Bag	7.5in (19.1cm)	9in (22.9cm)	1500ml	12mm	1/pkg, 5/cs	SB979

Customer Service: 1-800-448-6506
International Sales: 1-315-797-8375
Fax: 1-800-438-3051

ConMed.com
info@mail.conmed.com